ISTRUZIONI PER L'USO
INSTRUCTIONS FOR USE
MODE D'EMPLOI
BEDIENUNGSANLEITUNG
FOLLETO DE INSTRUCCIONES
HASZNÁLATI ÚTMUTATÓ

MASSAGGIATORE ANTICELLULITE
ANTI-CELLULITE MASSAGER
APPAREIL DE MASSAGE ANTICELLULITE
ANTI-CELLULITE MASSAGEGERÄT
MASAJEADOR ANTICELULITIS
ANTICELLULÓZ MASSÁZSGÉP

m240

ATTENZIONE

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ATTENTION

READ THIS BOOKLET
CAREFULLY.
IT CONTAINS IMPORTANT
INSTRUCTIONS FOR A SAFE
INSTALLATION, USE AND
MAINTENANCE.

IMPORTANT INSTRUCTIONS
TO BE KEPT FOR
FUTURE REFERENCE

WARNINGS

> GB

Keep handy the guarantee, receipt of payment and these instructions for future use and reference.

After removing the packaging, check the appliance for any signs of damage.

In case of visible damage, do not use it and contact a qualified technician.

Keep all parts of the packaging out of the reach of children. The appliance is not a toy: it is an electrical device and must be treated with the necessary caution.

Before connecting the appliance, make sure that the power supply voltage matches the indication on the data plate.

If the plug and socket are not compatible, replace the socket with a suitable type, applying to skilled people.

Do not use adapters or extension cables which do not meet current safety standards or that exceed the current carrying capacity limits.

When you don't use the appliance disconnect it from the power supply and assure that the selector is switched off (position 0).

Do not tug the power cable or the appliance itself to remove the plug from the socket.

The use of any electrical appliance requires to respect basic rules, namely:

- Never touch the appliance with wet or humid hands.

- This appliance must not be used by individuals (including children) with diminished physical, sensory or mental abilities or who lack the experience or knowledge on how to operate it unless they are supervised or have been instructed on the use of the appliance by someone who is responsible for their safety. Children must be supervised to make sure that they do not play with the appliance.

- Do not leave the appliance exposed to weathering (rain, sun).

- Do not subject it to rough handling.

The power cable should be extended to its full length to avoid overheating.

In case of a breakdown or malfunctioning of the appliance, disconnect it at once. Do not attempt to repair it, but take it to a service centre approved by the manufacturer.

If you decide not to use this appliance any longer, disable it by cutting the power cable - of course, after disconnecting it from the socket.

If the power cable gets damaged, replace it only in a service centre approved by the manufacturer.

For safety reasons, the appliance should never be opened.

WARNINGS

> GB

If you suffer from medical problems involving your feet or legs, such as phlebitis, varicose veins, skin excoriations or some other problems, consult your doctor for frequency or applications, as well as possible contraindications.

Expectant mothers, people who suffer from diabetes and people who carry pacemakers must consult their doctor before using this appliance.


Use the massager not over than 15 minutes.

Do not massage the same part of the body for more than 5 minutes. The appliance is not homologated for medical purposes.

TECHNICAL INFORMATIONS

> GB

- 1 - Functioning switch.
- 2 - massage speed regulating knob
- 3 - Accessories.
- 4 - Accessory cover in cloth
- 5 - Infrared.

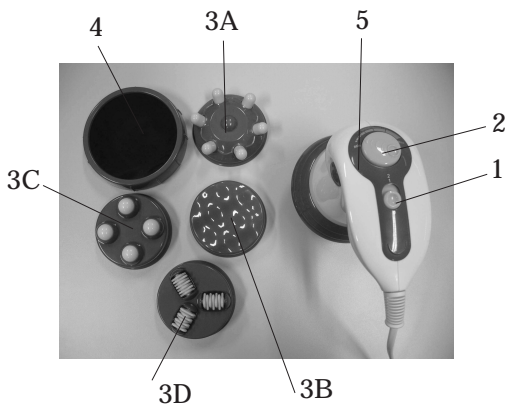This product has been manufactured in conformity to current 2006/95/EC and 2004/108/EC.

Technical data indicated on the appliance.

INSTRUCTIONS FOR USE > GB

Choose the accessory you wish to use for the massage, fit it to the hexagonal base making sure that the appliance is turned off and disconnected from the power outlet.

Fit the protecting cloth cap on the accessory 3B only and screw it onto the base with a clock-wise rotation.

Check to see that it is correctly closed.

The accessories provided are two (3A and 3B) to relieve muscular stress and two (3C and 3D) to reduce the "orange peel" appearance of the skin.

The infrared function helps the circulation of the blood during the massage.

Connect the plug into the power outlet and position the switch on 1 to activate the massage function alone, on 2 to add the infrared function.

Adjust the massaging speed using the specific knob.

Do not massage the same part of the body for more than 5 minutes and do not run the appliance for more than 15 minutes at a time.

Before beginning a new session of massages wait for the appliance to cool down at least 15 minutes.

MAINTENANCE > GB

This appliance does not require special maintenance. Clean it regularly. Before cleaning, make sure the appliance is switched off. Use a slightly damp cloth to clean the exterior of the appliance, while the accessories should be rinsed well and left to dry completely before being put away.

Use a soft, damp cloth to clean the appliance surface.

Do not use abrasive products or chemical solvents.

DO NOT IMMERSE THE APPLIANCE IN WATER.

IMPORTANT INFORMATION FOR CORRECT DISPOSAL OF THE PRODUCT IN ACCORDANCE WITH EC DIRECTIVE 2002/96/EC.

At the end of its working life, the product must not be disposed of as urban waste.

It must be taken to a special local authority differentiated waste collection centre or to a dealer providing this service.

Disposing of a household appliance separately avoids possible negative consequences for the environment and health deriving from inappropriate disposal and enables the constituent materials to be recovered to obtain significant savings in energy and resources. As a reminder of the need to dispose of household appliances separately, the product is marked with a crossed-out wheeled dustbin

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Terms and Conditions

The guarantee is valid for 24 months.

This guarantee applies only if it has been duly filled in and is submitted with the receipt showing the date of purchase.

The guarantee covers the replacement or repair of parts making up the appliance which were faulty at source due to manufacturing faults. After the guarantee has expired, the appliance will be repaired against payment.

The manufacturer declines any responsibility for damage to persons, animals or property due to misuse of the appliance and failure to observe the directions contained in the instructions.

Limits

All rights under this guarantee and any responsibility on our part will be voided if the appliance has been:

- mishandled by unauthorized persons
- improperly used, stored or transported.

The guarantee does not cover damage to the outward appearance or any other that does not prevent regular operation.

If any faults should be found despite the care taken in selecting the materials and in creating the product, or if any information or advice are required, please contact your local dealer.

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Prodotto importato da Ardes S.p.A.